



Single-piece hydrophilic acrylic  
Lenstec Softec HDM IOL

- » A thinner central optic which permits use of a CARTM cartridge allowing for a sub-2.2mm incision size.
- » Designed to be the most accurate IOL with patented bi-aspheric design with square edge technology
- » Tighter manufacturing diopter tolerance within +/- 0.11 D





## Micro Incision System Quarter Diopter • Bi-Aspheric

### Softec HDM Technical Specifications†

Optic Size	5.75 mm
Optic Type	Bi-aspheric
Length	12.00 mm
Haptic Style	Modified C
Angulation	0 Degrees
Positioning Holes	0
Construction	1 Piece
Optic Material	Acrylic Hydrophilic

†The A-Constant and ACD values printed are estimates only. It is recommended that the surgeon determine his or her own values based on their individual clinical experience.

### Constants (Optical Interferometry):\*

Immersion	A = 117.80
SRK/T	A = 117.80
Holladay I	sf = 1.11
Hoffer Q	pACD = 4.85

\*i.e. Using IOL Master, LENSTAR Optical Biometers  
\*Data Source: Early optimization (n=55) by N. Desai, T. Harvey.

### Diopter Steps


Whole	+5.00 to +36.00
Half	+10.50 to +29.50
Quarter	+15.00 to +25.00

## Softec HDM - Another Lenstec IOL designed to address spherical aberration and defocus.


**Softec HDM New Gold Standard:** 0.11 D Tolerance - 3x More Precise! Optical Prescription Selection and Tolerance Example: eye requires 24.25 D Prescription to achieve optimal vision.

Industry Standard IOL			Softec HDM IOL		
Lens	Allowed	Max. Variance for a standard IOL in an eye that requires a 24.25D lens is 0.65D. (smaller number is	Lens	Allowed	Max. Variance for a Softec HDM IOL in an eye that requires a 24.25D lens is 0.11D. (smaller number is
24.00	±0.4	)	24.00	±0.11	)
24.50	±0.4		24.25	±0.11	
			24.50	±0.11	
			24.75	±0.11	

**The Softec HD, Softec HDO, and Softec HDM are the only IOLs designed to address both Spherical Aberration and Defocus.** Defocus is a more significant aberration than Spherical Aberration.<sup>1</sup>

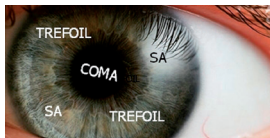

**Bi-Aspheric Equal Conic Zero aberration IOL.** Softec HDM addresses the issue of spherical aberration inherent in conventional monofocal spherical IOLs by adjusting the optic with a patented design on both the anterior and posterior surfaces.

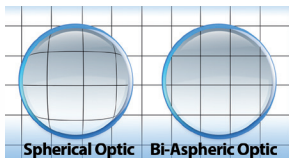

Studies have shown that Aspheric IOLs provide patients with significant optical benefits over traditional spherical surface IOLs.<sup>1,3,5</sup>

Softec HDM "Zero" Aberration:

- » Equal Bi-Aspheric
- » Less sensitive to decentration or tilt<sup>3,4</sup>
- » Ideal for all corneal profiles<sup>3</sup>
- » Enhanced depth of vision<sup>2</sup>

**Significant Outcomes.** The Softec HDM has been shown to help achieve refractive outcomes closer to intended<sup>6</sup>, significantly improve depth of field and decrease critical print size required for reading<sup>2</sup>, compared to a standard monofocal IOL.

### Proven quality - FDA approved

Lenstec is one of eight companies in the world approved by FDA (Food and Drug Administration) to sell intraocular lenses in the U.S. market. All products have CE certificate, are approved by BSI (British Standards Institute) and are ISO quality system certified.

Stability of the biomaterial from which the intraocular Lenstec lenses are made, is proven by a long term study and confirmed by millions of implanted lenses worldwide.

1. Thibos L, Hong X, Bradley A, Chang X. Statistical variation of aberration structure and image quality in a normal population of healthy eyes. *J Opt Soc Am A Opt Image Sci Vis* 2002; 19(12): 2329-48. » 2. Craig JP, Shah S, Wolffsohn JS. Clinical evaluation of the Softec HD aberration-free aspheric intraocular lens. *Clin Experiment Ophthalmol* 2011; 39(3): 281-3. » 3. Sarver E. Theoretical optical performance of an equal conic intraocular lens and comparison to spherical and aspheric IOLs. AAO Presentation 2005 » 4. Johansson B1, Sundelin S, Wikberg-Matsson A, et al. Visual and optical performance of the Akreos Adapt Advanced Optics and Tecnis Z9000 intraocular lenses: Swedish multicenter study. *J Cataract Refract Surg* 2007; 33(9): 1565-72. » 5. Nanavaty MA1, Spalton DJ, Boyce J, et al. Wavefront aberrations, depth of focus, and contrast sensitivity with aspheric and spherical intraocular lenses: fellow eye study. *J Cataract Refract Surg* 2009; 35(4): 663-71. » 6. Zudans JV, Desai NR, Trattler WB. Comparison of prediction error: labeled versus unlabeled intraocular

**World Headquarters**  
Lenstec, Inc.  
1765 Commerce Ave. N.  
St. Petersburg, FL 33716  
USA

Tel: 727-571-2272  
Fax: 727-571-1792  
Email: lenstec@lenstec.com

**Lenstec Barbados**  
Lenstec Barbados Inc.  
Airport Commercial Centre  
Pilgrim Road, Christ Church  
BARBADOS BB17092

Tel: 246-420-6795  
Fax: 246-420-6797  
Email: lenstecbarbados@lenstec.com

**Lenstec UK**  
Lenstec Barbados Inc.  
Lenstec House  
Unit 8, Mariner Court  
Calder Park, Wakefield  
West Yorkshire WF4 3FL  
England

Tel: +44 (0)1924 382 678  
Fax: +44 (0)1924 850 454  
Email: lenstecuk@lenstec.com



www.lenstec.com